



HAM10PD Atom A50B

10 Litre B & S Class Autoclave



The Atom HAM10PD is an ultra-fast compact autoclave that uses dry vacuum. It adopts **imported high-quality parts** including a **German** vacuum pump and an **Italian** water pump. The **system** ensures optimal drying, while the **patented independent steam generator** guarantees thorough sterilisation. The autoclave features **intelligent load detection** for energy efficiency, a **large storage capacity** for up to **30,000** records, and comprehensive safety measures including **automatic water quality monitoring** and emergency stop functions. The HAM10PD delivers reliable, **high-quality sterilisation** in a **user-friendly design**.



Fast Cycle Time From 12 Minutes



9 Programs: Cycle Parameters Can Be Customised To Meet Various Needs



User Friendly Design: Intuitive design combined with a 3.5" LCD display.



Large Storage Capacity: Can Store Up To 30,000 Sterilisation Records



Intelligent Load Detection: Calculates & Uses The Least Amount of Resources

Additional Features


- 3x Pre-Vacuum
- Pre-Heat Feature
- External Water Tank
- Intelligent Load Detection
- Intelligent Error Detection
- Embedded Control System
- Water Quality Monitor
- Built-In Large Memory
- 3x Trays 300x190x23mm
- Software Upgrade with USB
- Custom Program
- Emergency Functions
- Power Off Protection
- Preset Timer
- Auto Safety Door
- Easy Clean Jog-Dial Knob

High Quality Imported Parts

The series boasts imported high-quality parts, including a German Thomas vacuum pump (that can achieve a vacuum level of -0.90bar) and Italian solenoid valves and water pump, provide superior sterilization and maximum preservation for your instruments.

Optional External Printer

An optional Hatmed cycle record external printer can be easily connected at any time. No service technician needed.

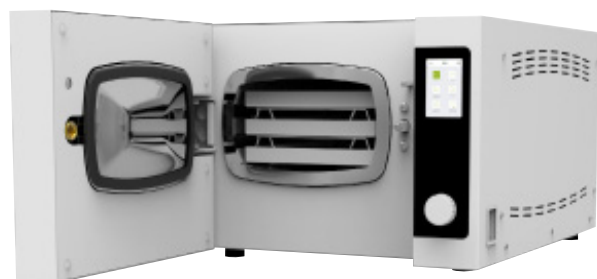

Cycle Durations (1kg Load & Including Drying Time)

| GEN134 | Quick B | Quick S | user1 | user2 |
|--------|---------|---------|-----------|-----------|
| 23 min | 14 min | 11 min | as needed | as needed |

Visit Our Website

