




HAM23PD Q70B

23 Litre B & S Class Autoclave

The HAM23PD is a compact air-cooling autoclave that uses dry vacuum. This is the air-cooling basic model, and the HAM23PDX is the pro model.

The HAM23PD adopts **imported high-quality parts** including a **German** vacuum pump and an **Italian** water pump. The **air-cooling system** ensures optimal drying, while the **patented independent steam generator** guarantees thorough sterilisation. The autoclave features **intelligent load detection** for energy efficiency, a **large storage capacity** for up to **30,000** records, and comprehensive safety measures including **automatic water quality monitoring** and emergency stop functions. The HAM23PD delivers reliable, **high-quality sterilisation** in a **user-friendly design**.

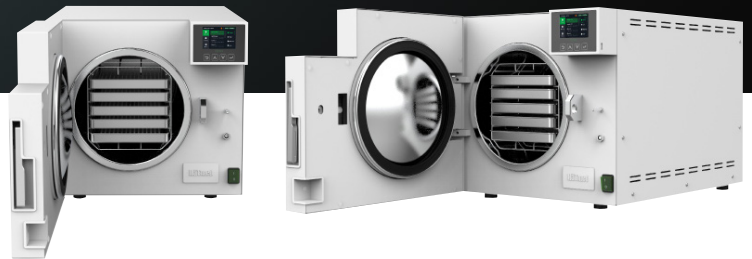

Fast Cycle Time From **26 Minutes**



9 Programs: Cycle Parameters Can Be Customised To Meet Various Needs



User Friendly Design: Intuitive design combined with a 2.4" LCD display.



Large Storage Capacity: Can Store Up To 30,000 Sterilisation Records

Additional Features

- 3x Pre-Vacuum
- Vacuum Drying
- Extra Dryness
- Intelligent Load Detection
- Embedded Control System
- Water Quality Monitor
- Intelligent Error Detection
- 5x Trays 395x190x23mm
- Built-In Large Memory
- Software Upgrade with USB
- Manual Lock
- Emergency Functions
- Power Off Protection
- Easy Clean Jog-Dial Knob

Optional External Printer

An optional Hatmed cycle record external printer can be easily connected at any time. No service technician needed.



High Quality Imported Parts

The series boasts imported high-quality parts, including a German Thomas vacuum pump (that can achieve a vacuum level of -0.90bar) and Italian solenoid valves and water pump, provide superior sterilization and maximum preservation for your instruments.

Cycle Durations (1kg Load & Including Drying Time)

GENI34	Quick B	Quick S	user1	user2
41 min	32 min	23 min	as needed	as needed

Visit Our Website

