



Ganj HAMG46TPD Q82BA

46 Litre B & S Class Autoclave



The 46L Ganj Autoclave is a large chamber steriliser that uses dry vacuum. The HAMG46TPD features imported high-quality parts including a **German** vacuum pump and an **Italian** water pump. Its **efficient condensing system** ensures optimal drying, while the **patented independent steam generator** guarantees thorough sterilisation.

The autoclave also features intelligent **load detection** for energy efficiency, a **large storage capacity** for up to **30,000** records, and comprehensive safety measures including **automatic water quality monitoring** and emergency stop functions. The HAMG46TPD delivers reliable, **high-quality sterilisation** in a **user-friendly design**.



Fast Cycle Time From **30 Minutes**



9 Programs: Cycle Parameters Can Be Customised To Meet Various Needs



5.6" Coloured Touchscreen: Elevates Usability with a User-Friendly Design.



Large Storage Capacity: Can Store Up To 30,000 Sterilisation Records



High Quality Imported Parts: Ensures Superior Sterilisation

Additional Features

- 3x Pre-Vacuum
- Pre-Heat Feature
- Vacuum Drying
- Intelligent Load Detection
- Embedded Control System
- Water Quality Monitor
- Intelligent Error Detection
- Built-In Large Memory
- Steel Wire Basket and Shelf
- External Water Tank
- Custom Program
- Preset Timer
- Software Upgrade with USB
- Manual Lock
- Emergency Functions
- Power Off Protection


High Quality Imported Parts

The series boasts imported high-quality parts, including a German Thomas vacuum pump (that can achieve a vacuum level of -0.90bar) and Italian solenoid valves and water pump, provide superior sterilization and maximum preservation for your instruments.

Cycle Durations (1kg Load & Including Drying Time)

| GENI34 | Quick B | Quick S | user1 | user2 |
|--------|---------|---------|-----------|-----------|
| 55 min | 36 min | 22 min | as needed | as needed |

Visit Our Website

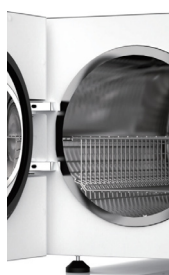

Optional Trays

The HAMG46TPD comes default with a flat steel wire shelf and a U-shaped steel wire basket. If you require more space to lay items on, an optional set of 5 trays can be purchased separately.

Optional External Printer

An optional Hatmed cycle record external printer can be easily connected at any time. No service technician needed.

